

# PHILIPS

# sonicare

## Power Flosser

## Gingival Health

in-vivo study

# A randomized, parallel study to compare Philips Sonicare Power Flosser and interdental brush on reducing gingival bleeding and plaque

Amini P, Imtiaz U, Li J, Ward M, Nelson M, Foster J

Silverstone Research Group, Las Vegas, NV, USA, 2020

Data on file

### Objective

The objective of this study was to compare the effects of manual toothbrushing plus two different interdental cleaning modalities, compared to manual toothbrushing alone, on gingivitis, gingival bleeding and plaque levels following a home use period of four weeks, and to assess the safety of the test products.

### Methodology

This was a three-arm, IRB-approved, single-blind, randomized, parallel clinical trial conducted in a generally healthy population of adults in the United States. Eligible subjects were aged 18–65 years, non-smokers and with no or infrequent interdental cleaning habits (less than once per week). Subjects were to have a minimum plaque score of at least 0.6, as scored by the Rustogi Modified Navy Plaque Index (RMNPI) following a three- to six-hour plaque accumulation period, as well as a gingival bleeding score of 1.0 on at least 50 sites, as scored by the Gingival Bleeding Index (GBI).

Eligible subjects were randomized to one of three treatment groups: manual toothbrushing alone (MTB), manual toothbrush plus TePe® interdental brush (IDB), or manual toothbrush plus Philips Sonicare Power Flosser with Quad Stream nozzle (PFQS). The Colgate Classic adult manual toothbrush and a standard dentifrice were dispensed to all groups and were used twice daily. Those randomized to IDB or PFQS performed adjunctive interdental cleaning once daily. The use of any other oral hygiene product or procedure was not permitted during the study period. Subjects were provided a diary to record compliance to the assigned regimen, as well as to document any safety events. Subjects returned to the clinic at Week 2 for an interim assessment of plaque and gingival health status, and at Week 4 for final evaluation, product return and dismissal from study.

### Results

There were 207 subjects screened, with 189 subjects enrolled and randomized to one of the three treatment groups: MTB (n=64), MTB+IDB (n=62), and MTB+PFQS (n=63). Of the 189 randomized subjects, 186 (98.4%) completed the study. There were 119 (63%) females and 70 (37%) males. The mean age of all randomized subjects was 37.2 (SD=12.5) years.

The overall distribution of race was white: 102 (54%), Black or African heritage: 42 (22.2%), Asian: 22 (11.6%), Native Hawaiian or other Pacific Islander: 13 (6.9%), and other: 10 (5.3%). For ethnicity, 44 (23.3%) of subjects randomized were Hispanic or Latino. There were no

significant differences in the treatment groups for demographic and baseline characteristics, including for plaque per RMNPI (p-value = 0.7351), gingival bleeding per GBI (p-value = 0.6738) and gingival inflammation per Modified Gingival Index, MGI, (p-value = 0.0952).

### Efficacy

#### Gingival Bleeding Index: LS Mean (SE) Baseline and Percent Reduction from Baseline Values

	MTB	MTB + IDB	MTB + PFQS
Baseline	0.50 (0.02)	0.51 (0.02)	0.49 (0.02)
<b>Percent Reduction from Baseline</b>			
Week 2	13.23% (2.25)	26.83% (2.31)	58.01% (2.23) **
Week 4	7.48% (2.19)	27.73% (2.21)	64.06% (2.18) **

\*Pairwise comparisons were statistically significant versus MTB treatment group, p-value < 0.0001

^Pairwise comparisons were statistically significant versus MTB + IDB treatment group, p-value <0.0001

#### Modified Gingival Index: LS Mean (SE) Baseline and Percent Reduction from Baseline Values

	MTB	MTB + IDB	MTB + PFQS
Baseline	2.73 (0.02)	2.77 (0.02)	2.71 (0.02)
<b>Percent Reduction from Baseline</b>			
Week 2	3.35% (0.60)	9.65% (0.62)	13.64% (0.60) **
Week 4	2.80% (0.63)	11.32% (0.64)	20.87% (0.63) **

\*Pairwise comparisons were statistically significant versus MTB treatment group, p-value < 0.0001

^Pairwise comparisons were statistically significant versus MTB + IDB treatment group, p-value <0.0001

#### Rustogi Modified Plaque Index: LS Mean (SE) Baseline and Percent Reduction from Baseline Values

	MTB	MTB + IDB	MTB + PFQS
Baseline	0.66 (0.00)	0.67 (0.01)	0.66 (0.00)
<b>Percent Reduction from Baseline</b>			
Week 2	5.28% (0.83)	10.99% (0.86)	20.42% (0.83) **
Week 4	6.86% (0.93)	14.47% (0.93)	28.78% (0.92) **

\*Pairwise comparisons were statistically significant versus MTB treatment group, p-value < 0.0001

^Pairwise comparisons were statistically significant versus MTB + IDB treatment group, p-value <0.0001

### Safety


There were 21 reported adverse events among 19 subjects, 14 of which were possibly related to the study, 2 in the MTB group, 10 in the MTB+IDB group, and 2 in the MTB+PFQS group. All events were characterized as mild in severity.

### Conclusions

All three treatment groups exhibited improvements in GBI, MGI, and RMNPI compared to baseline after two and four weeks of at-home use. The use of the Philips Sonicare Power Flosser, as an adjunct to manual toothbrushing, demonstrated the largest magnitude of reduction in gingival inflammation, bleeding and dental plaque, followed by manual toothbrush plus interdental brush, and manual toothbrush alone, at both two weeks and four weeks. Comparisons between products yielded statistically significant differences between the Philips Sonicare Power Flosser, the interdental brush, and manual toothbrush alone, for all measures, at both timepoints. All study products were safe for use.


**Figure 1:**  
**Gingival Bleeding Index**

LS Mean Percent Gingival Bleeding Reduction from Baseline to Week 2 and Week 4

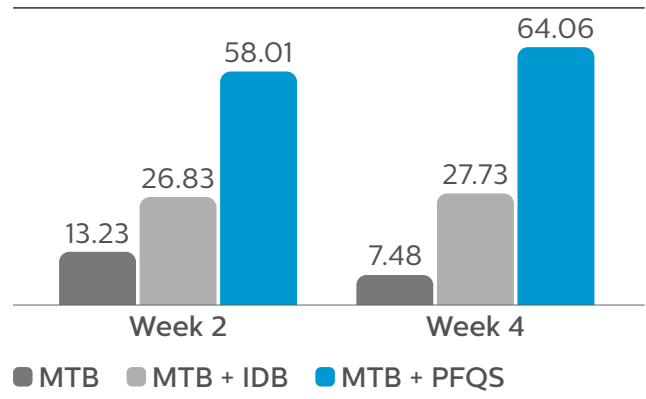

**Figure 2:**  
**Modified Gingival Index**

LS Mean Percent Gingivitis Reduction from Baseline to Week 2 and Week 4

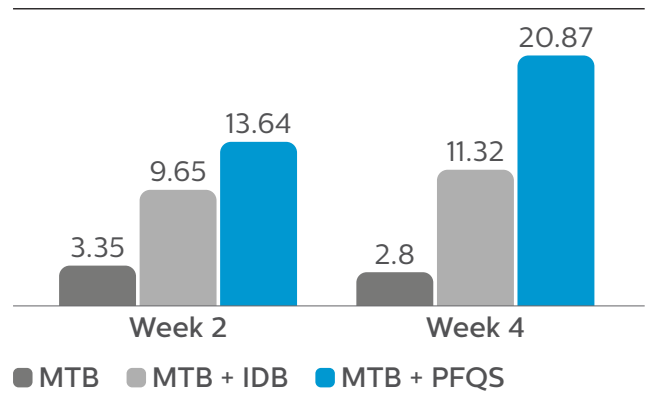

**Figure 3:**  
**Rustogi Modified Navy Plaque Index**

LS Mean Percent Plaque Reduction from Baseline to Week 2 and Week 4

