Turbo-Elite
Laser Atherectomy Catheter

Instructions for Use
For OTW and RX Catheter Models

Spectranetics

Turbo-Elite

Disposable Laser Optical Fiber Catheter
OTW and RX Catheters

Instructions for Use

Instructions for Use - Sections by Language

Pg	Language		Instructions For Use
3	English	English	Instructions For Use
15	Chinese	中文	使用说明

Table of Contents

1. PRODUCT DESCRIPTION	3
2. INDICATIONS	4
3. CONTRAINDICATIONS	5
4. WARNINGS	5
5. PRECAUTIONS.....	5
6. POTENTIAL ADVERSE EVENTS.....	6
7. CLINICAL RESEARCH	6
8. INDIVIDUALIZED THERAPY	8
9. OPERATOR MANUAL.....	8
10. HOW TO STORE.....	9
10.1 Sterilization.....	9
10.2 Check before Use	9
11. COMPATIBILITY	9
12. INSTRUCTIONS FOR USE.....	9
13. WARRANTY INFORMATION.....	13

US Federal law restricts this device to sale or on the order of a physician who has received appropriate training.

1. PRODUCT DESCRIPTION

Spectranetics Turbo-Elite disposable laser optical fiber catheter is a percutaneous intravascular appliance which contains a fiber bundle that forms a circle by surrounding the guidewire lumen plate.

For the Turbo-Elite Over-The-Wire (OTW) catheter, a luer connector adapter is installed at the proximal end of effective length so as to use the laser catheter installed with a guidewire of correct size (0.014”, 0.018”, and 0.035”). See the figure below.

For Turbo-Elite Rapid Exchange (RX) catheter, the guidewire lumen coaxial with the fiber bundle is located at the last 9 cm of the distal end and can come into direct contact with patient. See the figure below.

Action Mechanism of the Turbo-Elite catheter

The multifiber laser catheter transmits ultraviolet energy from the Spectranetics CVX-300 excimer laser system to the arterial obstruction. The ultraviolet energy is sent to the laser catheter tip and performs light resection of Fibrous, calcified and atherosclerotic lesions, thus probing the vessel with lesions (during light resection, high energy photons break molecular bonds at the cellular level without causing thermal damage to surrounding soft tissues). Spectranetics laser catheter is provided with a special lubricant coating and pass through the artery more easily.

Glossary of Special Terms

Retrograde mode = Opposite to the blood flow direction

Anterograde mode = Blood flow direction

Baseline angiography = Angiography record

Contralateral approach = Access to the artery using the crossover method

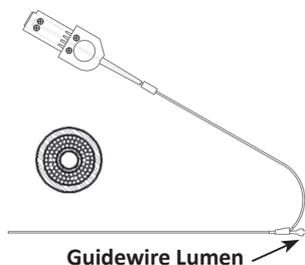
Figure 1: Turbo-Elite disposable laser optical fiber catheter (OTW type)

Table 1.1 Turbo-Elite disposable laser optical fiber catheter model

Device description	Model	Maximum guidewire compatibility (in.)	Maximum tip outer-diameter (in.)	Maximum shaft diameter (in.)	Working length (cm)	Sheath compatibility (Fr.)
OTW catheter specifications						
0.9 mm	410-152	0.014	0.038	0.047	150	4
1.4 mm	414-151	0.014	0.055	0.056	150	5
1.7 mm	417-152	0.018	0.068	0.069	150	5
2.0 mm	420-006	0.018	0.080	0.081	150	6
2.3 mm	423-001	0.018	0.091	0.091	125	7
2.5 mm	425-011	0.018	0.101	0.102	112	8
2.3 mm	423-135-02	0.035	0.091	0.091	125	7

Figure 2: Turbo-Elite disposable laser optical fiber catheter (RX type)

Table 1.2 Turbo-Elite disposable laser optical fiber catheter models

Device description	Model	Maximum guidewire compatibility (in.)	Maximum tip outer-diameter (in.)	Maximum shaft diameter (in.)	Working length (cm)	Sheath compatibility (Fr.)
RX catheter specifications						
0.9 mm	410-154	0.014	0.038	0.049	150	4
1.4 mm	414-159	0.014	0.057	0.062	150	5
1.7 mm	417-156	0.014	0.069	0.072	150	6
2.0 mm	420-159	0.014	0.080	0.084	150	7

2. INDICATIONS

The device is indicated for use with the Spectranetics CVX-300 excimer laser system and applies to the recanalization therapy of severe lower extremity arterial stenosis and occlusive lesions. The vascular reference diameter should be equal to or greater than 2 mm.

Note: When the guidewire is pushed step by step successfully, alleviation of severe limb ischemia cannot be guaranteed.

3. CONTRAINDICATIONS

No known contraindications.

4. WARNINGS

US Federal law restricts this device to sale or on the order of a physician who has received appropriate training.

Requirements for the CVX-300 excimer laser system software of Spectranetics Turbo-Elite disposable laser optical fiber catheter:

Software	Maximum repetition rate of catheter
V3.8XX	80 Hz
V3.7XX	40 Hz

When the laser catheter is in human body, use fluorescence microscopy equipped with X-ray photographic apparatus to perform observation so as to get HD images.

The laser catheter should not be operated in the presence of contrast. Prior to use, flush all residual contrast media from the introducer sheath or guide catheter and in-line connectors. For the 1.4-2.0 RX design, failure to do so may result in damage or dislodgement of the tip. Saline must be infused throughout the entire lasing process.

Fluoroscopic monitoring should be performed during use of the laser atherectomy catheter with special attention to identification of possible detachment of the marker band from the catheter tip. If detachment of the marker band is observed by fluoroscopy, it is recommended to leave the guidewire in place to aid in retrieval of the marker band, which can be attempted using balloon catheter inflation or snaring with the guidewire in place. If these techniques are unsuccessful, then stenting can be attempted to secure the marker band in place to avoid distal embolization. Please ensure you report all occurrences of this issue to Philips via a complaint.

5. PRECAUTIONS

This catheter has been sterilized using ethylene oxide and is a sterile product. This device is designed as a **disposable product**, so it cannot be re-sterilized or reused.

Do not re-sterilize or re-use this device, because incorrect re-treatment may hurt performance of the device and increase the risk of cross-infection.

Using this disposable device repeatedly may lead to serious casualties of patients and invalidates the manufacturer’s warranty.

Store the device at a cool and dry place. Avoid direct sunlight and high temperature heat source (exceeding 60 degrees Celsius or 140 degrees Fahrenheit).

The product sterility can be ensured only when the package is not opened and not damaged. Before use, check the sterile package carefully, and make sure that the seal is not torn. If the package integrity has been damaged, do not use this catheter. If the expiry date on the package label has passed, do not use this catheter product.

Before use, carefully check whether all the devices to be used during operation have any defects. Do not use the device if it has been damaged.

After use, dispose of all the devices according to the special requirements for hospital waste and hidden bio-hazardous materials.

Before operating the CVX-300 excimer laser system, please read through the sections related to warnings and responsibilities to be assumed in the operator manual carefully. These sections elaborate the notes, precautions and warnings to be observed to ensure the safe operation of the system.

 Spectranetics	Turbo-Elite Disposable Laser Optical Fiber Catheter OTW and RX Catheters	Instructions for Use
--	---	-----------------------------

English / English

During the procedure, provide the patient with the correct anticoagulant and vasodilator therapy according to the PTA protocol of the medical institution.

6. POTENTIAL ADVERSE EVENTS

Spectranetics CVX-300 excimer laser system may lead to the following complications:

Events observed in the clinical research (see section 7)

Surgery complications	Serious adverse events	Complications during hospitalization
<ul style="list-style-type: none"> • Spasm • Aortic dissection • Thrombus • Distal embolization • Perforation • Others 	<ul style="list-style-type: none"> • Death • Reintervention • ALI • Major amputation • Bypass grafting • Hematoma operation 	<ul style="list-style-type: none"> • Reocclusion • Renal failure • Pseudoaneurysm • Bleeding

Potential adverse events not observed in the clinical research (see section 7)

- | | |
|--|--|
| <ul style="list-style-type: none"> • Nerve injury • Formation of AV fistula • Endarterectomy • Infection | <ul style="list-style-type: none"> • Apoplexia • Myocardial infarction • Arrhythmia |
|--|--|

Since peripheral excimer laser recanalization is adopted, it is currently learned that long-term negative effects will not be caused to the arterial vascular wall.

7. CLINICAL RESEARCH

The data in this IFU are obtained from one subset of patients who come from three sources, and all have Critical Limb Ischemia (CLI). Therefore, they received continuous treatment, and operation is not suitable for them:

- LACI stage 2 – Patients for prospective IDE registration at 14 sites in America and Germany from 2001 to 2002 formed a subset. This subset included 26 lower limbs with lesions (25 patients) for treatment at 7 sites in America and Germany, and the stepping laser probing technology was adopted for all. Among these cases, 13 used the stepping technology in the very beginning, namely, the guidewire didn't get round the infarction as an attempt in the beginning.
- LACI Belgium – The 51 patients for prospective registration at 6 sites in Belgium formed a subset. This subset included 9 lower limbs with lesions (9 patients) for treatment at 3 sites in Belgium, and the stepping laser probing technology was adopted for all.
- Louisiana case series – A professional physician team continued data compilation at the Cardiovascular Institute of the South (CIS) in the middle of Louisiana, and selected 62 cases to form a subset. This subset included 12 lower limbs with lesions (12 patients), and the stepping laser probing technology was adopted for all.

 Spectranetics	Turbo-Elite	Instructions for Use
	Disposable Laser Optical Fiber Catheter OTW and RX Catheters	

English / English

Table 7.1. Procedure information

Arterial lesion location (n = 205)	
SFA	138 (67%)
Posterior knee	23 (11%)
Inferior genicular artery	42 (20%)
Results of arteriography (n= 47 cases of lower limbs)	
Lesions of each lower limb	4.4
Average damage length	73.4 ± 7.3 (mm)
Linear flow used for foot fixation	37 (79%)
Implanted stents	28 (60%)
Total number of successful penetration*	37 (79%)
Penetration success after guidewire insertion	24/34 (71%)
Cases with successful penetration counted from the beginning	13/13 (100%)
Operation success**	34 (72%)

Note: Totally 47 lower limbs with lesions of 46 patients were treated. The calculation of all the percentages was based on 47 cases.

* Traditional guidewires were inserted into 24 lower limbs. After an initio calculation was adopted for 13 lower limbs, the penetration success data of stepping cases was graded.

** Operation success: 50% terminal residual stenosis.

Table 7.2. Complications (n = 47 cases of lower limbs)

Surgical Complications	
Spasm	1 (2%)
Aortic dissection	4 (9%)
Thrombus	1 (2%)
Distal embolization	3 (6%)
Perforation	3 (6%)
Others	5 (11%)
Complications during hospitalization	
Reocclusion	1 (2%)
Aneurysmal hematoma	1 (2%)
Renal failure	1 (2%)
Bleeding	1 (2%)
Infection	0 (0%)
Others	0 (0%)

Note: Totally 47 lower limbs with lesions of 46 patients were treated. The calculation of all the percentages was based on 47 cases.

Table 7.3. Accumulated Serious Adverse Events (SAEs) in the follow-up period of 6 months (n = 47 cases of lower limbs)

Death	3 (6%)
MI or apoplexia	0 (0%)
Interventional therapy	6 (13%)
ALI	1 (2%)
Large resection	2 (4%)
Bypass surgery	2 (4%)
Endarterectomy	0 (0%)
Hematoma operation	2 (4%)
Total	16 (34%)

**Spectranetics****Turbo-Elite**
Disposable Laser Optical Fiber Catheter
OTW and RX Catheters

Instructions for Use

English / English

Note: Totally 47 lower limbs with lesions of 46 patients were treated. The calculation of all the percentages was based on 47 cases.

MI = Myocardial Infarction ALI = Acute Limb Ischemia

Table 7.4. Intentional analysis results (n = 47)

Penetration success	37 (79%)
Operation success	34 (72%)
Limb salvage	40 (85%)
Death, for any cause	3 (6%)
Any SAE	16 (34%)

NOTE: Totally 47 lower limbs with lesions of 46 patients were treated. The calculation of all the percentages was based on 47 cases.

8. INDIVIDUALIZED THERAPY

Before each patient uses the Turbo-Elite device, all of the above risks and benefits must be balanced carefully.

When traditional guidewire penetration fails in the beginning, the CLiRpath device can be considered because:

- The stump of circular or eccentric occlusion made the guidewire inclined to the vascular subintimal pathway.
- The guidewire was offset repeatedly and entered the large collateral branch flush with the stump of occlusion.
- Calcification blocked the guidewire path formed in the blocking cavity.

Select patients and apply clinical technologies according to instructions provided in section 2 “Indications” and section 9 “Operator Manual”.

9. OPERATOR MANUAL

The device operated in this document can be operated within the energy range of the CVX-300 excimer laser system:

Table 9.1 Energy Parameters

Device description	Model	Energy density	Repetition frequency	Laser on / off time
OTW Catheter				
0.9 mm	410-152	30-80	25-80*	Continuous On*
1.4 mm	414-151	30-60	25-80*	Continuous On*
1.7 mm	417-152	30-60	25-80*	Continuous On*
2.0 mm	420-006	30-60	25-80*	Continuous On*
2.3 mm	423-001	30-60	25-80*	Continuous On*
2.5 mm	425-011	30-45	25-80*	Continuous On*
2.3 mm	423-135-02	30-60	25-80*	Continuous On*
RX Catheter				
0.9 mm	410-154	30-80	25-80*	Continuous On*
1.4 mm	414-159	30-60	25-80*	Continuous On*
1.7 mm	417-156	30-60	25-80*	Continuous On*
2.0 mm	420-159	30-60	25-80*	Continuous On*

Recommended calibration settings: 45 energy density, 25 Hz.

* The maximum repetition frequency of software V3.812 is 80 Hz. The maximum repetition frequency of software V3.712 is 40 Hz.

10. HOW TO STORE

10.1 Sterilization

This device is for single use only. Do not re-sterilize or re-use this device.

The Spectranetics disposable laser optical fiber catheter is a sterilized product. Its sterilization performance can be guaranteed only when the package is not opened and not damaged.

10.2 Check before Use

Before use, check the sterile package, and make sure that the seal is not opened. Carefully check all the devices including the catheter to be used during the operation to see whether any defect exists. Check the laser catheter for bending, twisting or other damages. Do not use the device once it is damaged.

11. COMPATIBILITY

- The Spectranetics disposable laser optical fiber catheter is specially used with the Spectranetics CVX-300 excimer laser system.
- Do not use it with other laser systems.
- Guidewire Compatibility
- See the catheter specification table in section 1.

12. INSTRUCTIONS FOR USE

Installation procedure

The laser catheter package does not include the following additional materials, but some or all of these materials may be used during the procedure (all these materials are for single use only – do not re-sterilize or reuse them):

- Guiding sheath and femoral artery guiding catheter with the correct size and configuration, used to select peripheral artery and facilitate use of large laser catheters.
- Tuohy-Borst “Y”-adapter or hemostatic valve
- Sterile saline solution
- Standard comparison media.
- 0.014” and 0.018” guidewires

Only the radiologists with experience of peripheral vascular intervention can use the CVX-300 excimer laser system, and they must have received the required training, including but not limited to:

1. Safe laser emission.
2. Evaluate a patient’s injuries based on the film and according to the use instructions.
3. Evaluate various situations of applying the laser ablation technology to embolism according to the use instructions.
4. Evaluate the laser operation performed using the CVX-300 excimer laser system.
5. Personally practice using the CVX-300 excimer laser system of proper model.
6. Good practice of representative Spectranetics products will help to handle the first three situations.
7. Organized training together with exercise by physicians or exercises using the Spectranetics products through personal support will ensure efficiency of Spectranetics products.

Using sterilization techniques to open the aseptic package. Take down the packaging wedge from the bracket, lift the laser catheter from the bracket gently, and support the black laser connector at the same time, which is also called proximal adapter or proximal interface. Note that the proximal end of laser catheter can only be connected to the CVX- 300 excimer laser system, but cannot come into the contact with patient.

**English / English**


Connect the proximal end of laser catheter to the CVX-300 excimer laser system, and fix the laser catheter to the laser system extension rod. Calibrate the laser catheter according to the instructions for use in the operation manual of CVX-300 excimer laser system.

1. Use the standard femoral artery puncture technique to insert a guiding sheath of 4Fr. To 9Fr. (depending on the largest interventional device used during treatment) into the common femoral artery, and perform contralateral approach in the antegrade or retrograde manner. Adopt the heparinized PTA method to perform intravenous heparinization.
2. Inject contrast agent through the guiding sheath or guiding catheter, and carry out the baseline angiography. Use a composite projection device to obtain images, and describe the anatomic variation and morphology of the diseased region.
3. Insert a 0.014", 0.018" or 0.035" guidewire into the peripheral obstruction through the guiding sheath or guiding catheter.
4. Select a proper laser catheter according to the dimension arrangement:

Table 12.1 Recommended dimensions

Catheter dimension	The closest vessel diameter
0.9 mm	≥1.4 mm
1.4 mm	≥2.1 mm
1.7 mm	≥2.6 mm
2.0 mm	≥3.0 mm
2.3 mm	≥3.5 mm
2.5 mm	≥3.8 mm

5. Effect hydration on the catheter housing to activate the hydrophilic coating. Or immerse the catheter in a basin containing suitable sterile solution, or use wet gauze to wipe the catheter.
6. Use 5-10 ml of heparinized saline or lactated Ringer's solution to rinse the guidewire lumen of the laser catheter.
7. Use the selected guidewire to drag the distal end of the Spectranetics laser catheter. Use a fluoroscope to conduct monitoring, and guide the laser catheter to the lesion. The radiopaque marker of laser catheter displays the corresponding position of the catheter and the lesion.

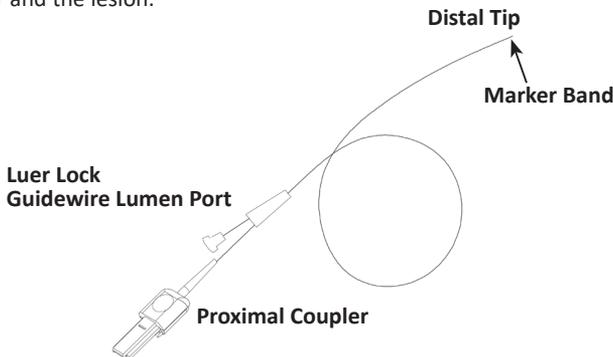**Figure 3 (not based on proportion)**

Note: During the use in human body, usually a fluoroscope is used to monitor moving of the laser catheter and position of the radiopaque tip position, similar to other intravascular interventional devices. The movement and advance rate of the distal end of catheter should be directly corresponding to advance rate of the upper shaft of catheter.

English / English

If corresponding movement is not apparent, reassess the lesion morphology, the laser energy being applied and the status of support equipment prior to continued treatment. If the catheter movement is not significant, operate more carefully, without excessively delivering the laser energy.

8. Inject contrast agent through the guiding sheath or guiding catheter, and use the fluoroscope to confirm the laser catheter position.
9. After making sure that the laser catheter has come into contact with the target lesion, use normal saline or lactated Ringer's solution to:
 - a. Rinse the residual contrast agent from the guiding sheath or guiding catheter and the internal connector.
 - b. Before activating the CVX-300 excimer laser system, rinse the residual contrast agent from the laser emission point and adjacent vascular structure.
Warning: Do not activate laser in the presence of contrast.
 - c. Perform saline solution flushing and infusion operation by referring to the content related to the saline infusion method in the use instructions.
10. Once the foot switch is pressed during use of the Turbo-Elite catheter, the CVX-300 excimer laser system will release energy continuously. The operator can control the laser chain length. Generally, the recommended duration of continuous laser emission cannot exceed 20s.
11. **Stepping method of treating total occlusion**
 - a. Press the foot switch to activate the CVX-300 excimer laser system, and then slowly insert the laser catheter into 2-3 mm at the position of total occlusion, enabling the laser energy to resect the target substance. Release the foot switch to stop the CVX-300 excimer laser system.
 - b. Drag the guidewire beyond the distal end of laser catheter, further insert it into the occlusion by about several millimeters, and then activate the laser again according to the description in the above step 8.
 - c. Continue the operation in the stepping manner to make the guidewire and laser catheter move forward and activated (with the increment of millimeter), till the catheter reaches the last 3-5 mm of occlusion.
 - d. When passing through the last 3-5 mm of the occlusion, make the guidewire enter the distal vessel first, and then the activated OTW laser catheter.
 - e. Fix the guidewire in place, withdraw the laser catheter, use the guiding catheter to inject contrast agent, and monitor the lesion through the fluoroscope.
 - f. To remove more lesion, the additional laser transmittance test needs to be performed for the OTW catheter.
 - g. If the catheter meets an obstacle (e.g., calcium) when moving forward, release the foot switch immediately to interrupt operation of the CVX-300 excimer laser system and stop laser emission. To move forward the catheter continuously, increase the energy density and repetition frequency. To avoid heat accumulation, the catheter must be pushed during laser emission.
12. **Standard method for stenosis treatment**
 - a. Press the foot switch to activate the CVX-300 excimer laser system, and then make the laser catheter pass through the stenosis slowly at a speed smaller than 1 mm per second. Release the foot switch to stop the CVX-300 excimer laser system.
 - b. To remove more lesion, the additional laser transmittance test needs to be performed for the OTW catheter. If the catheter meets an obstacle (e.g., calcium) when moving forward, release the foot switch immediately to interrupt operation of the CVX-300 excimer laser system and stop laser emission. To move forward the catheter continuously, increase the energy density and repetition frequency. To avoid heat accumulation, the catheter must be pushed during laser emission.

**English / English**

13. To increase/reduce the energy density and pulse repetition frequency, the laser catheter does not need to be removed from the patient's body because the laser catheter has been calibrated beforehand. Refer to the operation manual of CVX-300 excimer laser system.
14. If necessary, the follow-up angiography and balloon angioplasty should be performed after the laser probing. A stent can be placed for cases such as acute contraction and severe perforation according to requirements.
15. Saline infusion method.

Note: This technique requires two operators. It is recommended that the main physician – operator should drag the catheter, and control the foot pedal of laser system. The assistant responsible for wiping should use the control syringe of saline infusion, and press the fluoroscope pedal (if applicable).

- a. Before the laser surgery, heat one bag of 500 ml 0.9% normal saline (NaCl) or lactated Ringer's solution to 37°C. Heparin or potassium does not need to be added to the normal saline/LR solution. Connect this bag of hot normal saline/LR solution to the sterile venous catheter, and connect one end of the catheter to the triplicate three-way connector.
- b. If applicable, use a conventional method to insert an appropriate guiding catheter with a "large lumen" into the opening of artery. It is advised to avoid side holes for the guiding catheter.
- c. Push the laser catheter under guidance of the fluoroscope, till it comes into contact with the lesion. If necessary, inject contrast agent to help to position the tip of laser guide. If contrast agent stays between the tip of laser guide and the lesion, the laser catheter needs to be slightly retracted (by 1 to 2 mm) so that normal saline/LR solution can implement anterograde flow and remove contrast agent when being used to flush the system. Before laser emission, however, make sure that the tip of laser guide has come into contact with the lesion.
- d. Discharge the remaining contrast agent from the control syringe, and place it back to the contrast agent bottle. Suck the normal saline/LR solution into the control syringe through a three-limb tube, and clear the contrast agent from the triplicate three-way connector.
- e. Remove the control syringe preliminarily installed on the three-limb tube, and replace it with a new control syringe of Luer lock with a capacity of 20 ml. Before connection, this new 20 ml controllable syringe should suck normal saline/LR solution to reduce the occurrence probability of bubbles. (Merit Medical and other manufacturers who produce 20 ml control syringes).
- f. Use at least 20-30ml of normal saline/LR solution (several syringes of saline/LR) to rinse the blood and contrast agent from the three-limb tube, connector pipe, Y type connector and the guiding sheath or guiding catheter. After the first rinsing process is completed, inject normal saline/LR solution into the 20 ml control syringe again.
- g. Use the fluoroscope to confirm that the tip of laser catheter has come into contact with the lesion (push to move the catheter when necessary). Do not inject contrast agent here.
- h. When the main operator indicates that he/she has got ready for activation of the laser system, the assistant responsible for wiping should turn the three-limb tube piston rapidly to eliminate its pressure, and inject 10 ml of normal saline/LR solution (within 1 to 2 seconds). This kind of bolus injection aims to transfer the blood or dilute it to the level of blood capillary, and restrict blood reflux of the laser ablation area.

**English / English**

- i. When the injection activity still continues after the initial 10 ml bolus injection is performed, the assistant responsible for wiping should slow down the injection to 2-3 ml/second. This proportion of normal saline/LR solution injection aims to transfer or dilute the antegrade blood flow entering the laser ablation area. Now, the assistant responsible for wiping should press the foot switch properly to activate the CVX-300 excimer laser system and start emitting the laser beam.
- j. The operator controls the laser chain length. Generally, the recommended duration of continuous laser emission cannot exceed 20 seconds. The normal saline/LR solution injection must last throughout the laser emission process.
- k. Stop injecting normal saline/LR solution when the laser chain ends. Turn the three-limb tube piston to make it feed back pressure, and suck 20 ml of normal saline/LR solution again using the control syringe, getting ready for the next time of laser beam emission.
- l. Before emitting each beam of laser chain subsequently, use normal saline/LR solution to perform bolus injection, and inject normal saline/LR solution continuously according to the description in steps h-k.
- m. If contrast agent is used to evaluate the therapeutic effect during laser therapy, reactivate the CVX-300 excimer laser system (repeat steps h-k before activating the excimer laser system).

Note: Normal saline/LR solution can be injected through the sheath (anterograde method) or laser catheter lumen (contralateral method) according to the used method, antero- or contralateral. When the contralateral method is used, it is advised to use a guidewire with a smaller diameter to inject normal saline/LR solution into the procedure site thoroughly.

13. WARRANTY INFORMATION

Limited warranty of the manufacturer

The manufacturer ensures that the Turbo-Elite disposable laser optical fiber catheter does not have any defect in materials and processes provided that it is used in the validity period and its package has not been opened or damaged before use. According to this warranty, the manufacturer's responsibility is limited to replacement of the defective Turbo-Elite disposable laser optical fiber catheter, or refunding according to the purchase price only. The manufacturer shall not be responsible for any accidental, special or consequential damages due to use of the Turbo-Elite disposable laser optical fiber catheter. This limited warranty shall become invalid if the Turbo-Elite disposable laser optical fiber catheter is damaged due to misuse, transformation, improper storage or operation, or other faults occur due to the failure of performing operations according to the use instructions. **This limited warranty can be in lieu of all the other explicit or implicit warranties, including any implied merchantability or fitness for a particular purpose.** Any individual or entity, including the agent or intermediary authorized by the manufacturer, shall not have the right to expand or extend this limited warranty, and any attempt made out of this purpose shall have no binding force to the manufacturer. This limited warranty covers the Turbo-Elite disposable laser optical fiber catheter only. The manufacturer's warranty information for this system can be found in the related documents of CVX-300 excimer laser system.

**Spectranetics****Turbo-Elite**
Disposable Laser Optical Fiber Catheter
OTW and RX Catheters

Instructions for Use

English / English

Series	Model	Guidewire compatibility (mm)	Outer diameter of tip (mm)	Maximum outer diameter of catheter (mm)	Transmission efficiency (%)	Effective length (cm)
OTW	410-152	0.356	≤0.965	≤1.19	≥2.94	150
	414-151	0.356	≤1.40	≤1.42	≥6.89	150
	417-152	0.457	≤1.73	≤1.75	≥11.56	150
	420-006	0.457	≤2.03	≤2.06	≥17.28	150
	423-001	0.457	≤2.31	≤2.31	≥27.78	120
	425-011	0.457	≤2.57	≤2.59	≥41.67	110
	423-135-02	0.89	≤2.31	≤2.31	≥20.28	120
RX	410-154	0.356	≤0.965	≤1.24	≥2.94	150
	414-159	0.356	≤1.50	≤1.57	≥6.50	150
	417-156	0.356	≤1.75	≤1.83	≥9.44	150
	420-159	0.356	≤2.03	≤2.13	≥16	150

[Sterilization Method] Ethylene oxide sterilization

[Validity Period] 2 years

[Manufacturer] Spectranetics Corporation

[Address of registrant / Manufacturer] 9965 Federal Drive Colorado Springs
Colorado 80921 USA

[Contact method of registrant / manufacturer] 1-800-231-0978

[Manufacturing address] 9965 Federal Drive Colorado Springs Colorado 80921 USA

[Agent / After sales service organization name] Philips (China) Investment Ltd. Co.

[Address of Agent / After sales service organization] Building A1, #718 Lingshi Road,
JingAn District, Shanghai

[Contact method of Agent / After sales service organization] 8008100038

[Production Date] See the Label.

[Product Technical requirement Serial Number] Imported Medical Instrument
20163012942

[Medical Device Registration Certificate No.] Imported Medical Instrument 20163012942

[Environmental Conditions] Storage and transportation temperature: not exceeding 60°C;
storage and transportation humidity: avoid moisture, keep dry

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©2025 Spectranetics Corporation
Manufactured by Spectranetics Corporation